


# HOW TO TEST YOUR DRUGS FOR *FENTANYL*

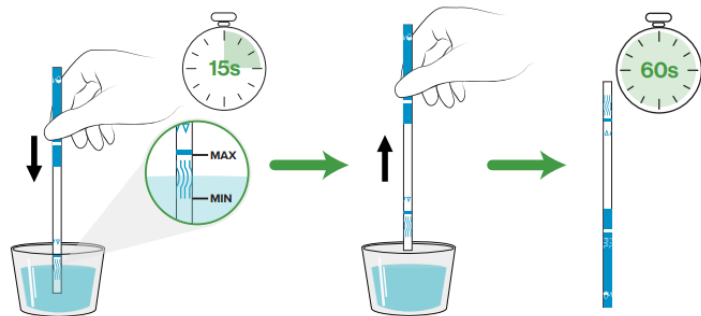

Strips cannot detect all types of fentanyl, such as carfentanil.

If using alone, call

**National Overdose Response Service**

1-888-688-NORS (6677)

For more information, visit  
**timiskaminghu.com**



Services de santé du  
**TIMISKAMING**  
Health Unit



1. Put your drugs (5-10 mg or a tiny pinch) into 5 ml (teaspoon) of room temperature water.
2. Stir until the drug is completely dissolved.
3. Dip the strip into the cup for **15 seconds**, no deeper than the wavy line.
4. Wait 1 minute for **red lines** to appear. Do not read after 10 minutes.
5. Check the test results. **One line means that the substance has fentanyl in it.**



**1 LINE:**

**POSITIVE** for fentanyl



**2 LINES:**

**NEGATIVE** for fentanyl

A negative test result does not mean the substance is safe. **Always use caution. Carry Naloxone.**